



Thermal Desorption Sorbent Tubes

For Active Sampling of ppb-level VOCs in Ambient or Indoor Air



KEY FEATURES

- High sensitivity
- High recovery
- Reusable

Thermal Desorption

SKC sorbent technology research produces superior materials for the collection of organic compounds for thermal desorption.

Thermal desorption techniques offer the advantage of greatly improved analytical sensitivity. Because a solvent is not used for extraction, the collected sample is not diluted and, in most cases, analytical recovery is so close to 100% that desorption efficiency corrections are not required.

SKC Sorbents Well Suited for Thermal Desorption

To be suitable for thermal desorption, sorbents must meet exacting specifications that include:

- Low contaminant background
- High thermal stability
- Sufficient adsorptive strength to retain components of interest, but release them quickly when heat is applied

SKC sorbent technology provides superior materials for the collection of organic compounds for thermal desorption.

Anasorb® 747 (20/40 mesh size and 980 m²/gm surface area)

This beaded active carbon has the capacity for organic vapors similar to petroleum-based and coconut-shell charcoals. It also can be used to sample nonpolar organic compounds.

Anasorb GCB1 (20/40 or 60/80 mesh size and 100 m²/gm surface area)


Anasorb Graphitized Carbon Pack Black 1 (GCB1) has a moderate surface area and is equivalent to Carbpac B (60/80). This sorbent has proven to be valuable for sampling compounds of intermediate to high volatility.

Anasorb GCB2 (20-40 or 60/80 mesh size and 10 m²/gm surface area)

This sorbent has a low surface area and is equivalent to Carbpac C (60/80). Anasorb GCB2 is useful in collecting semi-volatile organic compounds (SVOCs).

Carbpac X (40/60 mesh size and 240 m²/gm surface area)

This strong graphitized carbon sorbent collects C3 to C9 hydrocarbons. Carbpac X is hydrophobic.

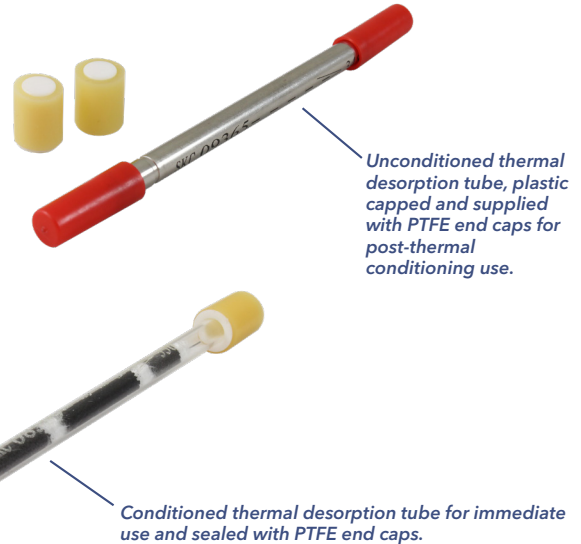

Carbosieve S-III (60/80 mesh size and 975 m²/gm surface area)

This carbon-based molecular sieve has properties similar to the obsolete Anasorb CMS and replaces it in SKC thermal desorption tubes. Carbosieve S-III has a high capacity/breakthrough volume for small molecules and is moderately hydrophilic. It can be used for the lower boiling compounds such as halogenated hydrocarbons and alcohols.

Chromosorb® 106 (60/80 mesh size and 700 to 800 m²/gm surface area)

A porous polymer sorbent, Chromosorb 106 is used to trap low-boiling hydrocarbons, benzene, labile compounds, and volatile oxygenated compounds. This hydrophobic sorbent is the least polar polymer in the Chromosorb family and is functionally similar to Anasorb 727.

Tenax® TA (35/60 or 60/80 mesh size and 30 to 35 m²/gm surface area)

Tenax TA is a traditional sorbent (porous polymer) for trapping medium to high boiling compounds; it is especially useful for low concentrations because of its low background. Tenax TA is hydrophobic and is suitable for use in EPA Method TO-17 or IP-1B and other thermal desorption applications.

Tenax GR (60/80 mesh size and 35 m²/gm surface area)

This Tenax/graphite composite extends the range of Tenax to lower boiling compounds. It retains the ability of Tenax to be cleaned easily to a very low background and is widely used in thermal desorption applications.



SCIENCE. SERVING PEOPLE.®

www.skinc.com | The Leader in Sampling Solutions and Expertise for OEHS Professionals

More information 

2 | Thermal Desorption Sorbent Tubes

Thermal Desorption Tube Selection for Perkin Elmer or Markes International Thermal Desorbers

SKC offers single and multiple-bed thermal desorption tubes that meet EPA Method TO-17 requirements. SKC thermal desorption tubes are available conditioned for immediate use and sealed with PTFE end caps. SKC also offers unconditioned thermal desorption tubes that are plastic capped and supplied with PTFE end caps for use after thermal conditioning by your laboratory as needed. Each SKC thermal desorption tube is marked with a permanent serial number. It is best practice to condition thermal desorption tubes before each use to ensure background levels are suitable for ppb-level measurements.

ORDERING

Perkin Elmer or Markes International Thermal Desorber Tubes

Available in **glass or stainless steel**, these tubes measure 1/4-inch OD x 3 1/2-inch length (6.35-mm OD x 88.9-mm length). Supplied as a single tube.

Method	Sorbent	SS Conditioned Cat. No.	SS Unconditioned Cat. No.	Glass Conditioned Cat. No.	Glass Unconditioned Cat. No.
ASTM D6196-15e1	Anasorb GCB1*	226-356	226-356-UP	–	–
ASTM D6196-15e1 MDHS 104†	Tenax TA	226-357	226-357-UP	226-360	226-360-UP
ASTM D6196-15e1 MDHS 104†	Chromosorb 106	226-358	226-358-UP	–	–
EPA TO-1, IP-1B	Tenax TA	226-340	226-340-UP	226-339	226-339-UP
EPA TO-2	Carbosieve S-III	226-341	226-341-UP	–	–
EPA TO-17	Anasorb GCB1*/Carbosieve S-III	226-349	226-349-UP	–	–
EPA TO-17, NIOSH 2549	Anasorb GCB2*/Anasorb GCB1*/Carbosieve S-III	226-350	226-350-UP	226-347	226-347-UP
EPA TO-17	Tenax GR/Anasorb GCB1*	226-348	226-348-UP	–	–
EPA TO-17	Anasorb GCB2*/Carbosieve S-III	–	–	226-346	226-346-UP
EPA TO-17	Tenax GR/Anasorb GCB2*	–	–	226-345	226-345-UP

*Anasorb GCB1 is equivalent to Carbopack B; Anasorb GCB2 is equivalent to Carbopack C.

†Stainless steel tubes only

Perkin Elmer or Markes International Thermal Desorber Tubes with Carbopack X

Select tubes with Carbopack X sorbent for superior recovery of very volatile compounds. Available in stainless steel, these tubes measure 1/4-inch OD x 3 1/2-inch length (6.35-mm OD x 88.9-mm length). Supplied as a single tube.

Method	Sorbent	Cat. No.	Qty.
MDHS104†	Carbopack X	Conditioned	226-363
		Unconditioned	226-363-UP

†Stainless steel tubes only

Screening Tubes (Dynatherm Thermal Desorber)

Available in **glass**, these tubes measure 0.24-inch OD x 4.5-inch length (6-mm OD x 115-mm length). Supplied as a single tube.

Method	Sorbent	Cat. No.	Qty.
EPA IP 1B, NIOSH 2549	GCB2*/GCB1*/Carbosieve S-III	Conditioned	226-330
		Unconditioned	226-330-UP

*Anasorb GCB1 is equivalent to Carbopack B; Anasorb GCB2 is equivalent to Carbopack C.

Accessories

	Cat. No.	Qty.
Glass Transport Tubes, for 3.5-inch (88.9-mm) length tubes	226-300	5
Swagelok® Fittings, for 1/4-inch (6.35-mm) OD tubes	P50291	2
PTFE Ferrules	P30121	2

**For additional details, visit
Thermal Desorption Tubes on our website.**